




Pressure Care Cushions

Pressure care cushions are designed to alleviate the pain and discomfort often associated with spending prolonged periods of time in fixed positions. From high-density foam to dry floatation air cushions these products offer a range of solutions to suit specific needs.

NOVIS COMFORT CUSHION

Ideal to improve comfort of users in manual wheelchairs. A super-soft memory foam upper layer is supported by a high density foam base. A polyurethane 2-way stretch top cover provides waterproof and vapour permeable protection. The cover can be easily unzipped and removed for cleaning. W410mm x L410mm x H75mm

Product Code: *Nov-comf cush*

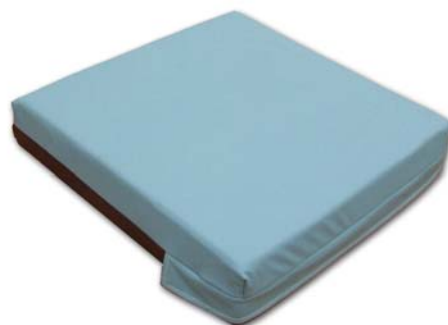

NOVIS PREMIUM CUSHION

A super-soft memory foam upper layer is shaped and sized for additional comfort and wheelchair use. A polyurethane 2-way stretch top cover and non slip base provides waterproof and vapour permeable protection. The cover can be easily unzipped and removed for cleaning. W430mm x L410mm x H105/65mm

Product Code: *Nov-prem wh cushion*



SYNERGY SIMPLICITY

The Synergy Simplicity low-profile cushion offers great comfort plus the health benefits of moderate pressure management and positioning. It features high density polyurethane foam construction for long lasting comfort, wrap-and-lift design for pressure management and proper positioning, coccyx and sacral relief and a 4-way stretch, fluid resistant, fabric cover. The cover can be easily unzipped and removed for cleaning.


Product Code: *P-CUSASMB+Size ?* DVA Code: *AE05 225603*



SYNERGY SOLUTION

For the active client, the Solution has a contoured, wrap-and-lift base designed with high density, high resiliency polyurethane foam that redistributes weight away from areas vulnerable to skin breakdown. This, along with twin-cell gel well, viscoelastic foam layer, 4-way stretch top fabric, pliable urethane fluid resistant inner liner and non slip bottom, offers pressure management, positioning and long lasting comfort. The cover can be easily unzipped and removed for cleaning.

Product Code: *P-CUSASMB+Size ?* DVA Code: *AE05 225606*

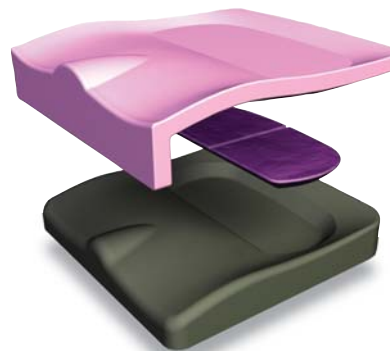




Pressure Care Cushions Continued..

SYNERGY SPECTRUM

The Spectrum takes the wrap-and-lift base design to the next level by providing an extra-deep well and allowing the provider to choose the medium that offers the best pressure management and comfort for the individual client. Choices available are: the lightweight twin-cell gel or the viscoelastic foam insert.

This system provides comfort as well as excellent pressure management and positioning. The cover can be easily unzipped and removed for cleaning.


Product Code: P-CUSASMB+Size ?

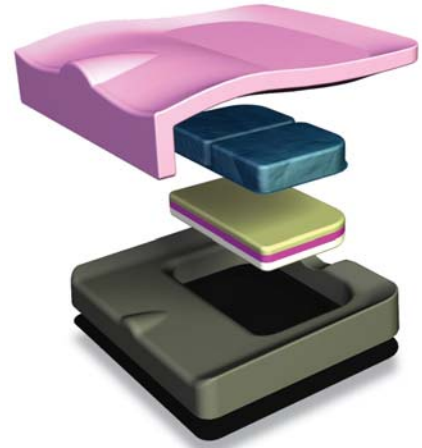

SYNERGY STRUCTURE

The Structure expands even further on the wrap-and-lift base design by combining a choice of pressure management inserts - the lightweight twin-cell gel or the viscoelastic foam - with modular positioning components including additional wedges and abductors. This allows us to deliver client-specific customisation for maximum pressure management, unparalleled positioning and incredible comfort. The cover can be easily unzipped and removed for cleaning.

Product Code: P-CUSASMB+Size ?

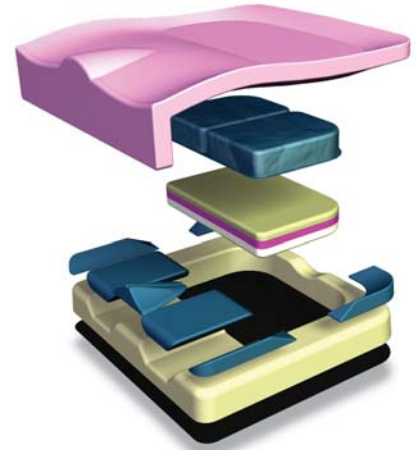

NEXUS SPIRIT CUSHION

The Nexus Spirit cushion combines the stability of a contoured foam base with the therapeutic care for which Roho Dry Floatation cushions are known. The contoured foam base is designed to allow increased stability for transfers and provides exceptional positioning of the pelvis and lower extremities for enhanced sitting posture.

The Roho Cushion, positioned to address critical pressure points, provides you the highest level of skin protection from the original Roho Dry Floatation Technology.

Product Code: SD-NEXUSSPIRIT+Size ? DVA Code: AE04 220442



HARMONY CUSHION

The Roho Harmony cushion combines a lightweight and comfortable Roho air flotation cushion with a pre-contoured foam base and is encased within a moisture resistant cover.

Particularly suited to users requiring stability, ease of transfer and postural support.

Product Code: SD-HARMONY+Size ?





Pressure Care Cushions Continued..

QUADTRO SELECT CUSHION

The Quadtro Select cushion sets the standard in the industry for overall performance in wheelchair seating. As the name suggests the cushion is divided into 4 compartments, each of which can be set for optimal comfort and support. Every user will benefit from a customised fit through the simple push of a knob. The revolutionary Isoflo Memory Control offers shape fitting capabilities while the user is seated, allowing quick and easy, on-demand adjustment to maximise function.

With the built-in stability and simplicity of the Quadtro Select cushion, there is no need to sacrifice maximum skin protection to get stability, positioning or convenience. Available in an extensive range of sizes.

Product Code: SD-QS+Size ? DVA Code: AE04 220444



JAY 2

The Jay 2 Cushion is the industry leading adjustable fluid cushion. Lightweight and low maintenance, the cushion is designed for users at high risk of skin breakdown. The pre-contoured foam base and 2" of 3D Jay Flow fluid ensure correct postural stability and comfort. Custom cuts and modifications are available to meet specific client needs. Available in an extensive range of sizes.

Product Code: J-J2+Size? DVA Code: AE04 21+Size?

JAY 2 DEEP CONTOUR

The Deep Contour cushion provides even greater comfort and support for users at extremely high risk of skin breakdown. The pre-contoured foam base and 3" of deep immersions 3D Jay Flow fluid tri pad ensure correct postural stability and comfort. Custom cuts and modifications are available to meet specific client needs. Available in an extensive range of sizes.

Product Code: J-J2D+Size? DVA Code: AE04 24+Size?



INVACARE FLOTECH

Highly adaptable gel cushion with a contoured foam base. Dual gel sacs cover the entire top surface of the cushion and re-distribute pressure away from bony prominences. An additional layer of gel is incorporated under the ischial and sacral area to further enhance the pressure reducing qualities of this cushion, making it suitable for Very High Risk users.

The cushions can be further enhanced with the Flo-Tech Modular System which further promotes active correction and posture management. Velcro pads on the cushion and the components ensure the wedges remain fixed in place even after extensive use. All modular components can be customised and are made from high-density CMHR foam.

Product Code: MT-MSS+Size? DVA Code: AE05 226103

